


Operating Manual
About Prospera

Prospera develops high quality, innovative, affordable massage products—from electronic pulse massagers to jade massagers—for active people. Prospera products are available through leading department stores, consumer electronics stores, and major online retailers. Learn more about Prospera at www.prosperacorp.com.

Prospera Corporation

405 Boulder Court, Suite 500
Pleasanton, CA 94566
Tel: (925) 225 0888
www.prosperacorp.com

Introduction

Prospera's PL022 Electronic Slim Massager is a user friendly product. Prospera's Electronic Slim Massager tones muscles in your waist, arms, legs and abdomen with its powerful knock and high-speed vibration massage. You can regulate the intensity of the massage with four operating modes: constant vibrating, pulse, wave, and cycle of vibrating.

Your muscles will sigh with relief when you soothe them with the Electronic Slim Massager. The ultra-light unit comes with a remote control and can be easily adjusted with a Velcro strap, to bring maximum sensation just where you need it. Firm up muscles at home or in the office!

PL022 comes with a transformer. PL022 is a perfect personal massager for the home and office.

Features of PL022 Electronic Slim Massager:

1. Multiple speeds with different intensities
2. Four operating modes: constant vibrating, pulse, wave, and cycle of vibrating.
3. Heat function accelerate blood circulate.
4. Gentle and comfortable massage
5. Well designed controller for easy and simple operation
6. Easily adjusted with a strap, to bring maximum sensation just where you need it.
7. A perfect personal massager for the home and office

Content
Content

Safety Warnings	3
Set up	4
Operating the product	5-6
Trouble shooting	7
Consumer warranty	8
Positions	9
Warranty card	10
Prospera company contact	

Safety Warnings

When the Electronic Slim Massager is not in use, always unplug it from the transformer. This action conserves energy, prevents electric shock and damage to the product.

Warnings:

Please take these precautions to prevent electric shock, burns, fires, and other injuries.

1. Do not leave this product unattended when it is in use. Unplug the product from the power surge protector or cigarette lighter after usage.
2. Do not operate this product in an excessively humid environment, and avoid direct sunshine upon the product where it is set up.
3. Do not insert any sharp objects into the fabric of the product. These objects might include keys, pens, pencils, kitchen utensils, blades, knives, etc.
4. Do not use the product outdoors.
5. Do not use air refresher spray where the product is set up or in use.
6. Consult your doctor prior to using the product if you have heart, skin and other health issues.
7. This product is not intended for children to use. Keep close supervision if your child uses it.
8. This transformer is recommended to be plugged into a surge protector to avoid damage to its mechanism and electric components.

Attention:

1. Do not let babies and young children to play with the product, remote, or cord.
2. If you feel that the massage is too strong, stop program immediately, and disconnect the product from the transformer.
3. Average daily usage for massage product is about 15 to 30 minutes. Avoid watching television when using this product, as you may be unaware of how much time has passed. Do not use it while you are asleep.
4. This product contains electric components, wires, and other elements that cannot contact excessive moisture. For instance, the bathroom, shower room, steam room, and sauna are places where this product should not be used.
5. Do not sit on the massager to massage.

3

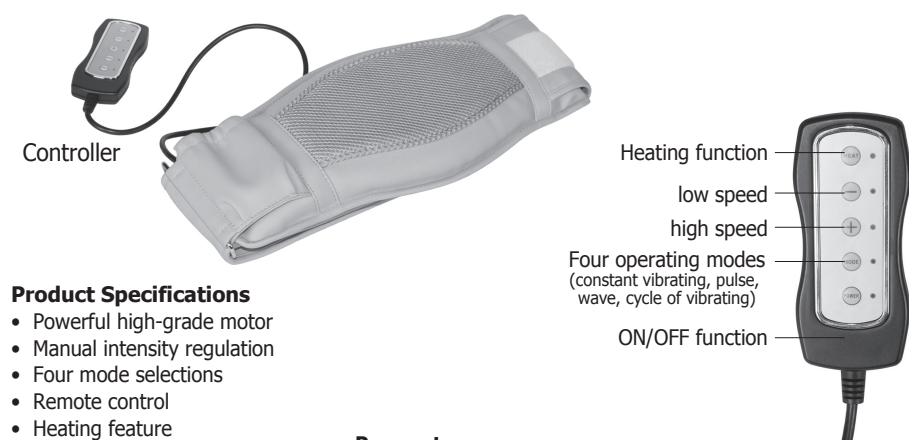
Set up

Find a good location in your living room, family room, study, or office to set up the Prospera's PL022 Electronic Slim Massager. Use a power surge protector for connecting the transformer that comes with the product.



Take the product out of retail box, place it at your waist, or shoulders, and insert the plug to the transformer and then to the surge protector. Turn the surge protector on.

4

Operating the product

Product Specifications

- Powerful high-grade motor
- Manual intensity regulation
- Four mode selections
- Remote control
- Heating feature
- Adjustable velcro strap
- Input 12V (adapter included)
- Weight: 2.25 lbs
- Measures: 17" x 7" x 4.5"
- UPC: 185277000219

Parameters:

Model: PL022
Name: Electronic slim massager
Input: DC12V, 1.2A
Power: 14.4W

Transformer:

Input: 120VAC, 60Hz
Output: DC12V, 1.2A

5

Operating the product
Function:

1. Multiple speeds with different intensities
2. Four operating modes: constant vibrating, pulse, wave, and cycle of vibrating.
3. Heat function accelerate blood circulate.
4. Gentle and comfortable massage.

Operating instruction:

Turn on the power, choose massage mode to start massage session. Press Heat if you prefer this feature. Manually regulate speed, "-" or "+" to get the most suitable speed and intensity for your massage session.

After you are done massaging, switch the Power off . Unplug the insert from the transformer.

Recommended time of usage: 15 to 30 minutes daily

6

Trouble shooting

Prospera's PL022 Electronic Slim Massager is a user friendly product. It is built to last. If you follow the operating manual, do not over use it, and follow the precautions, this product will last for a long time.

Should the following situations happen, please take recommended trouble shooting procedures:

1. If the fabric is torn, stop using it immediately.
2. If a wire is broken or the product does not work, unplug the product from the surge protector. Contact Prospera if it is under warranty.

7

Consumer Warranty

Prospera offers a 90-day standard warranty to the massage products. Labor and parts are included during warranty period. Consumers are required to provide model number and dated proof of purchase (sales receipt) when contacting Prospera for the warranty service.

Non-warranty service is provided on a "per incident" basis. The consumer will verify that the product has failed and provide information for servicing the unit. All applicable repairs, parts, shipping, handling, local tax and "per incident" fees will be charged for non-warranty repairs and support calls.

8

Positions


9

Warranty card

Warranty Registration Form

Product Name: _____
 Product Model: _____ Purchase Price: _____
 Date of Purchase: _____ Retailer: _____
 First Name: _____ Last Name: _____
 Address: _____
 City: _____ State: _____ Zip: _____
 Phone: _____ Email: _____
 Signature and date: _____

Please mail to:

Prospera Corporation Telephone: 925-225-0888
 405 Boulder Ct, Ste 500 Fax: 925-225-0660
 Pleasanton, CA 94566 Email: sales@prosperacorp.com



10