

AWS

OM-200

User Manual



OM-200 (440lb x 0.2lb)

OM-200

Thank you for purchasing the American Weigh OM-200 digital bathroom scale. Please read all operating instructions carefully before use. Scales are precision measuring devices and should always be handled with proper care. To ensure years of reliable service, keep these simple tips in mind:


- Do not exceed the scale's maximum capacity. Overloading your scale can permanently damage it!
- Avoid exposure to extreme heat or cold. Scales perform best at normal room temperature. If temperatures have changed dramatically, calibration and adjustment may be necessary.
- Allow your scale to warm up for 30-60 seconds before performing calibration and adjustment(if available).
- Store your scale in a clean, dry location. Dust, dirt, and moisture can accumulate on the weighing sensors and electronics causing inaccuracy or malfunction.
- Avoid using your scale in close proximity to strong radio wave emitting devices such as computers, cash registers, and cordless phones.
- Always weigh on a flat and level surface, free from vibrations and drafts. The corner of a room is usually the most stable.
- Gently apply all items to be weighed. Do not drop items onto the weighing platform.
- Avoid dropping your scale. The warranty does not cover damage due to rough treatment or overload.
- Check the batteries first if you are having any trouble with your scale. This simple step can remedy most scale issues.
- Do not disassemble your scale. This product contains no user serviceable parts.

I. Part List

1. Scale
2. CR2032 Battery (x1)


II. Key Description

Bottom of Scale

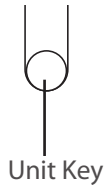

III. Battery Installation

1. Before using the scale for the first time, check the battery compartment and remove any plastic insulation strips or packaging that may be in place to prevent battery drainage during shipping/storage.
2. If you are replacing the batteries, make sure the positive and negative contacts are properly aligned.

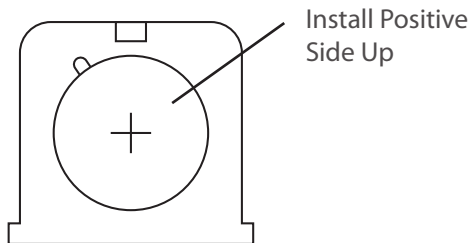

IV. Scale Operation

1. Place the scale on a flat, level surface. Do not weigh on carpet as it may affect the accuracy.
2. Step onto the platform. Once the weight has stabilized, it will be locked on the display for a few seconds.
3. The scale will turn off automatically after 10 seconds of non-use.

NOTES: To prevent slipping, do not use the scale with wet feet. Do not leave items on the scale when not in use. This scale has an Auto-On feature that allows you to step on the scale immediately without having to tap the scale first to zero it. If moving the scale from one place to another, allow the scale to turn on then off once before weighing. This will allow the zero-point to be re-established. Otherwise, the weight result may be inaccurate.

V. Changing the Unit of Measure

1. Press the **UNIT** key on the bottom of the scale to change the scale's displayed unit of measure.
 2. If the scale is turned off then back on, it will default to the last unit used.
- **OM-200 Units:** kg (Kilogram) ▶ lb (Pound) ▶ st:lb (UK Stone:Pound)

VI. Specifications

	OM-200
Max. Capacity	440lb / 200kg / 31st
Readability	0.2lb / 0.1kg / 1lb
Power	1 x CR2032
Dimensions	11.8 x 11.8 x 0.8"
Display Dimensions	2.1 x 3.6"
Warranty	10 Year Limited Warranty

American Weigh Scales, Inc.
Norcross, GA 30092, U.S.A

Made in China to our exact specifications and quality control.