

ENDURANCE®

— BUILT TO LAST! —



Overbed Table

Item #P2600



7 54756 62600 6

MADE IN TAIWAN


2011 ESSENTIAL MEDICAL SUPPLY, INC.®



Made for
ESSENTIAL®

MEDICAL SUPPLY, INC.
Orlando, FL 32822

PARTS

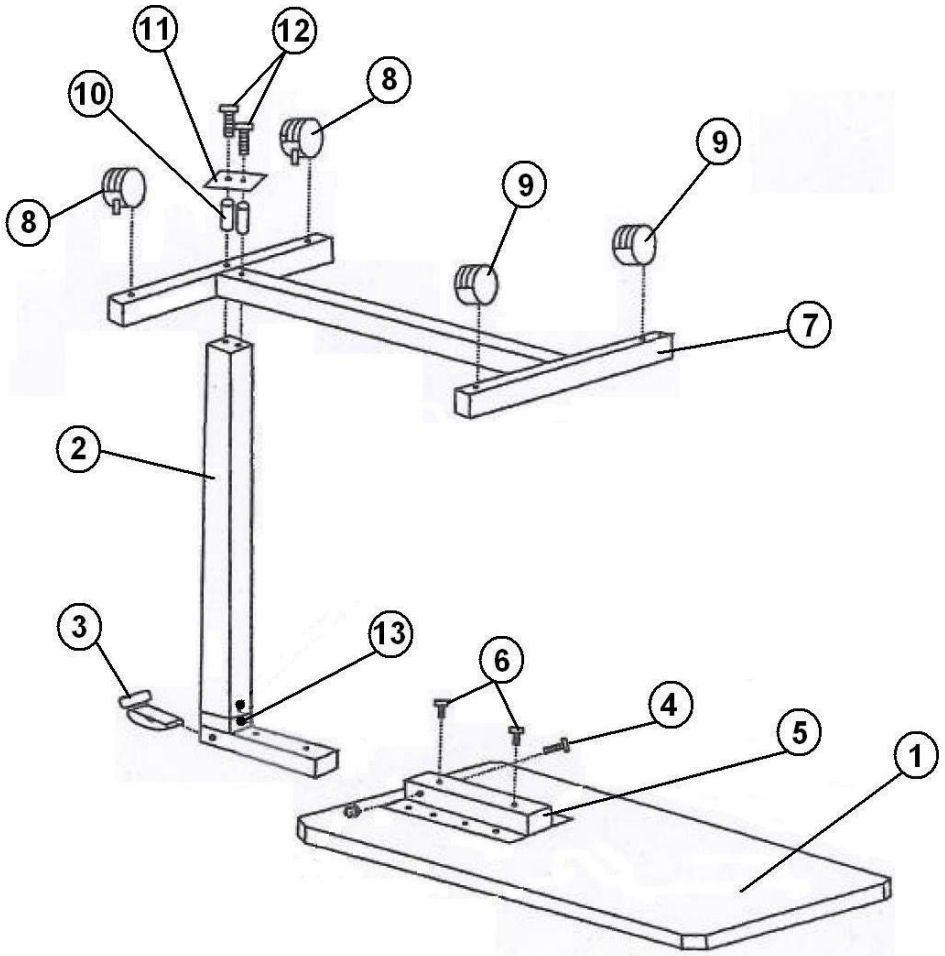

PARTS LIST

- | | | | |
|---|----|--------------------------------|----|
| 1. Table Top..... | x1 | 7. Table Base..... | x1 |
| 2. Column Assembly..... | x1 | 8. Caster with Lock..... | x2 |
| 3. Trip Lever..... | x1 | 9. Caster without Lock..... | x2 |
| 4. M6 x 1 x 50mm Long Hex Head Bolt & Hex Lock Nut..... | x1 | 10. Steel Bushing..... | x2 |
| 5. Table Mount Bracket..... | x1 | 11. Steel Plate..... | x1 |
| 6. M8 x 1¼ x 12mm Long Bolt..... | x2 | 12. M10 x 1½ x 50mm Bolt..... | x2 |
| | | 13. Black Flat Head Screw..... | x1 |

P2600 OVERBED TABLE

ASSEMBLY INSTRUCTIONS

1. Invert Table Top (1) and slip the short arm of the Column Assembly (2) into the Table Mount Bracket (5). Make sure the Trip Rod is fully inserted inside and secure with two M8 Bolts (6).
2. Install Trip Lever (3) into the Column with rounded side facing the Table Top (1). Align the holes with the Table Mount Bracket (5), insert M6 Bolt (4) and secure it with the Lock Nut (4).
3. Install Table base (7) onto the Column (2) and insert the two Steel Bushings (10) into holes in the table base (7), place the Steel Plate (11) and secure with M10 Bolts (12) using 14mm wrench.
4. Install two Casters (8) and (9) into the Table base (7), as shown in diagram.
5. Turn the Table assembly right side up and remove the Black Flat Head Screw (13) from the Column.
6. To adjust the Table Top (1) to desired height, squeeze the Trip Lever (3) slightly and pull or push down the Table Top. (1).

CAUTION

- SPRING LOADED COLUMN ASSEMBLY (2)
- Do NOT remove the Black Flat Head Screw (13) before the full Assembly is complete.
- The Casters with Locks (8) shall be in locked position before using the Table.