



Bedpan Mounting Bracket

Item 6820A100

Contents:

Bedpan Holder (1)

Spring Hook (2)

"L" Bracket (4)

7/8" Tap Screw (4)

1/2" Bolt (4)

Nylon Locking Nut (4)


Tools needed:

Phillips screwdriver

ASSEMBLY INSTRUCTIONS:

1. Remove seat assembly from the frame of the chair and set face down on a flat working surface.

Fig. 1

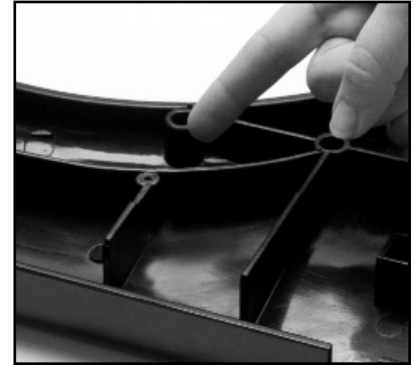

2. Locate four small screw holes (fig. 1) on seat assembly and install "L" brackets with a Phillips screwdriver and 7/8" tap screws (fig. 2).

Fig. 2



3. Drop one 1/2" bolt into one spring hook, ensuring that the threaded portion of the bolt is exposed out of the bottom of the spring (fig 3).

Fig. 3



4. Install spring hook onto the flat cross-section of the bedpan holder by inserting the threaded end through one hole and hand tighten using a nylon locking nut (fig. 4). Repeat step 3 and install on the opposite end.

Fig. 4



5. Insert the attached spring hooks through the exposed holes on the “L” brackets at the back of the seat assembly (fig. 5).

Fig. 5




6. Attach the rest of the bedpan holder to the two remaining “L” brackets by inserting a bolt through the hole on the bracket and the opening on the holder (fig. 6). Place a nylon locking nut over the threaded end of the bolt and screw the bolt into the nut until it locks.

Fig. 6



7. Re-install the seat assembly onto the frame of the chair, ensuring that the bedpan holder handle is placed below the back horizontal bar (fig. 7).


Fig. 7

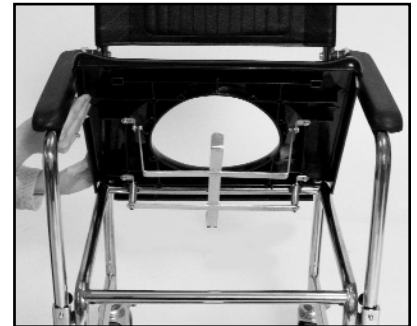

Bedpan Installation

Pull the bedpan holder handle down and insert bedpan into the holder assembly (fig. 8). Release the handle and the bedpan should rest snugly beneath the seat assembly.

Fig. 8

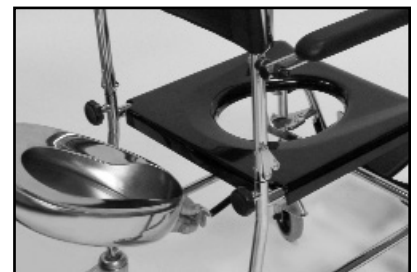

Graham-Field Health Products

2935 Northeast Parkway, Atlanta, Georgia 30360
telephone: 800-347-5678, 770-447-1609
fax: 800-726-0601, 678-291-3232
<http://www.grahamfield.com>
Made in China