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Escali
Measuring What's Next™

Model: BFBW200
User's Guide

Body Fat / Body Water Analyzing Scale

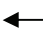
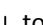
www.escali.com

Thank you for choosing the Escali Body Fat / Body Water analyzing scale. This precision state-of-the-art measuring instrument is most accurate when weighing on a flat hard surface. Avoid carpet. As the electronic sensors are sensitive, avoid dropping or jarring the scale and store where it will be protected from impact.

 **Warning!** Do not use this product if you have a pacemaker or other internal electronic medical device.

Initial Setup & Entering Personal Data Into Memory

1. Remove the insulating plastic strip from the batteries and ensure the +/- terminals are correctly placed as shown inside the battery compartment. This scale operates on 2 - CR2032 lithium batteries (included).
2. Select the desired unit mode with the Kg/St/Lb slide switch on the base of the scale.
3. Turn the scale on by pressing the  Set/Select function key.
4. Select user number (1 - 10) using  or , press  to confirm.
5. Select gender using  or , press  to confirm.
6. Select height using  or , press  to confirm.
7. Select age using  or , press  to confirm.
8. After completing programming, the LCD will return to zero. Your data is now stored into memory. Remember your memory location (1-10). You will need it each time you measure your body composition.
9. To program other users repeat steps 3 through 7.

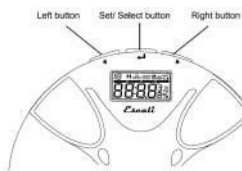
TO CHANGE BATTERIES

Ensure +/- terminals are correct way round.

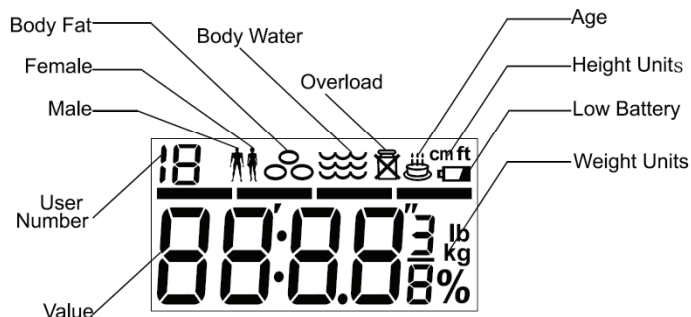
*Please remove the insulating strip before use.

FUNCTION KEYS

NOTE: KEYS ARE LOCATED ON THE SIDE OF THE SCALE.

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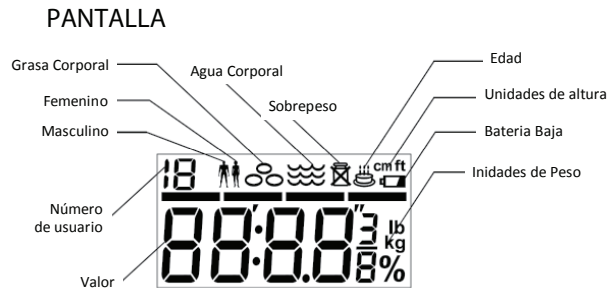
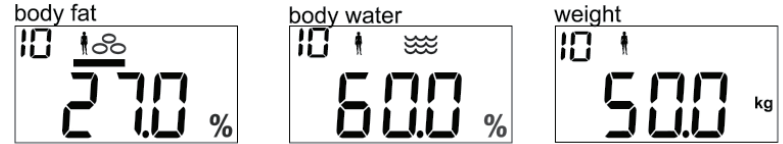
To Measure Body Fat & Water Content

1. Use your foot to press and hold the center part of the scale platform. You will see '1' is increasing, 1-10. When the scale displays your desired User Number, release your foot.
2. When the display reads '0.0', step onto the scale platform with **bare feet** and stand still.
3. Your body fat & water in % and weight will be shown on the display as follows one at a time.

4. The scale turns off automatically if no further operation.

To Measure Weight Only

1. Use foot to tap the platform scale.
2. When display shows '0.0', step onto the platform and stand still.
3. Your weight value is then shown on the display.
4. The scale turns off automatically if no further operation.

How does the BFBW200 work?

This scale uses BIA (Bioelectrical Impedance Analysis) to measure the body fat and body water. In this technique, a safe, low-level electrical signal is passed through the body. The scale measures the speed of the current which is affected by water, muscle, bone, and fat. Using this data, the scale estimates the percentage of body fat and water.

Using the Analysis

Hydration levels fluctuate throughout the day and are affected by exercise, eating and drinking. As a result, for maximum accuracy and repeatability, it is recommended that this function be used about the same time of day under the same conditions to establish a personal average. These readings are a close indication, and not an exact value. Please consult your doctor for more information.

Weight Loss

During weight loss, the percentage of body fat may not change as expected since the reading is a percentage of body weight and both will be changing at the same time. For example a 180lb person with 30% body fat has 54 lbs of fat (30% of 180 lb = 54 lb). If that person loses 15 pounds but still has 30% body fat, they have still lost 4 ½ pounds of fat (30% of 165 lb = 49.5 lb / 54 - 49.5 = 4.5). To minimize the loss of lean body tissue be sure to include exercise in your weight loss plan.

HEALTH TABLE

The Body Fat Ratio (%)

Age	Female			
	Underfat	Healthy	Overfat	Obese
10-39	0-21	22-33	34-39	40-80
40-59	0-23	24-34	35-40	41-80
60-79	0-24	25-36	37-42	43-80

Age	Male			
	Underfat	Healthy	Overfat	Obese
10-39	0-8	9-20	21-25	26-80
40-59	0-11	12-22	23-28	29-80
60-79	0-13	14-25	26-30	31-80

The Water Ratio (%)

Age	Female			Male		
	Low	Normal	High	Low	Normal	High
10-15	<57	57-67	>67	<58	58-72	>72
15-30	<47	47-57	>57	<53	53-67	>67
30-60	<42	42-52	>52	<47	47-61	>61
60-80	<37	37-47	>47	<42	42-56	>56

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Additional Information

1. Make sure your scale is at room temperature, and on a stable horizontal surface free of vibration.
2. This scale is not waterproof. Do not expose to moisture.
3. As this scale is a precision instrument, it must be handled with extreme care. Avoid rough treatment, shaking and vibration.
4. Do not disassemble this scale or the warranty will be void.
5. The display value may be affected by electromagnetic disturbances such as operating near a radio. If this occurs, move the device that is causing the disturbance and turn the scale off and restart.

Product Guarantee

Your Escali scale is guaranteed to be free of defects in materials and workmanship for a lifetime. For full details of this guarantee please visit our web site at: www.escali.com. Escali's liability is limited to two (2) times the cost of the product.

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