**Escali®**

Measuring What's Next™

Model: USTT200

User's Guide

# Track & Target Bathroom Scale

[www.escali.com](http://www.escali.com)

Thank you for choosing the Escali Track & Target Bathroom Scale. This precision state-of-the-art measuring instrument is most accurate when weighing on a flat hard surface. This scale will not function correctly on carpet. As the electronic sensors are sensitive, avoid dropping or jarring the scale and store where it will be protected from impact.

## Operating Instructions

1. Place four AAA batteries (total 4) into the battery compartment located on the base of the scale. Ensure the +/- terminals are correctly placed as indicated by the diagram within the battery compartment.
2. Select your desired mode of display (lb or kg) by pressing the button located on the base of the scale.
3. Press the Set Button. The word User will begin to flash. Using the Up and Down arrows, select the User number preferred. Press Set to confirm selection.
4. The word Target will begin to flash. Using the Up and Down buttons select the target weight. Press Set to confirm selection.
5. The word Weight will appear and display 0.0. The scale is now ready for use. Simply stand still on the weighing platform with your feet evenly spaced within 10 seconds. Once your weight has been determined, a message to Get Off the scale will be displayed. Once you have stepped off of the scale the progress toward the target weight will be displayed along with the change in weight since the last weighing.
6. If you do not step off of the scale when the words Get Off are displayed, the scale will begin to rotate through the 4 User Numbers. When the display shows your desired User Number step off the scale and your progress will be displayed.

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
7. To change the target weight simply repeat steps 3 5.
8. For future use simply stand directly onto the scale. As long as your weight has not changed more than 6.6 lbs (3 kg) the scale will select your user information automatically.

### Additional Information

1. Make sure your scale is at room temperature, and on a stable horizontal surface free of vibration.
2. This scale is not waterproof. Do not expose to moisture.
3. As this scale is a precision instrument, it must be handled with extreme care. Avoid rough treatment, shaking and vibration.
4. Do not disassemble this scale or the warranty will be void.
5. The display value may be affected by electromagnetic disturbances such as operating near a radio. If this occurs, move the device that is causing the disturbance and turn the scale off and restart.

### Product Guarantee

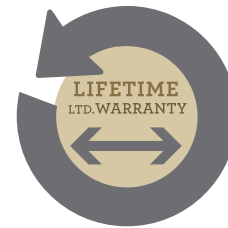Your Escali scale is guaranteed to be free of defects in materials and workmanship for a lifetime. For full details of this guarantee please visit our website at:

Escali's liability is limited to two (2) times the cost of the product.

[www.escali.com](http://www.escali.com)

**Escali®**

Measuring What's Next™

Modelo: Serie USTT200

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# Báscula Track & Target

[www.escali.com](http://www.escali.com)