Escali[®]
Measuring What's Next[™]

Model: XL200

User's Guide

**Bathroom Scale
Extra-Large Display**

www.escali.com

OPERATING INSTRUCTIONS

Thank you for choosing the Escali XL180, Extra Large Display scale.
This precision state-of-the-art measuring instrument is most accurate when weighing on a flat hard surface.

Avoid carpet. As the electronic sensors are sensitive, avoid dropping or jarring the scale and store where it will be protected from impact.

FEATURES & SPECIFICATION:

1. Capacity : 200 KG / 440 LB
2. Weight Graduation : 0.1 KG /
3. Operation Temperature : 0~40 °C
4. Power Supply: 2 x CR 2032 Battery
5. Weight Units : LB, KG
6. LCD V.A. : 95 mm x 75 mm
7. Auto Off In 6 Seconds If No Operation
8. Low Voltage Indication: "Lo"
9. Error Indication : "Err"

FULL DISPLAY

INITIAL SETUP :

1. Remove the insulating plastic strip from the batteries and ensure the +/- terminals are correctly placed as shown inside the battery compartment. This scale operates on 2-CR2032 lithium batteries (included).
2. Select desired unit mode (Kg or Lb) with slide switch on bottom of scale.

HOW TO WEIGH:

1. Use foot to tap on the middle of the weighing platform.
2. LCD shows a moving '8', then '0,0'.
3. Now step on the platform and stand still.
4. Your weight reading is shown in about 3 seconds then locked.
5. You can now step off the scale. The display will hold your weight reading for about 5 seconds
6. The scale then turns off automatically in about 8 seconds.

TURNING OFF YOUR SCALE:

The scale turns off automatically in about 8 seconds.

CHANGING THE WEIGHING MODE:

You can simply change the weighing mode on the back side of the scale.

DO NOT PLACE OVERWEIGHT ITEMS ON THE SCALE:


Weighing items over the maximum weight capacity will cause an overload 'Err' error to display. Immediately remove overweight items to avoid causing damage to your scale.

BATTERY REPLACEMENT:

Replace the batteries when the LCD display shows the low battery 'LO' symbol.

TO CHANGE BATTERIES

Ensure +/- terminals are correctly placed

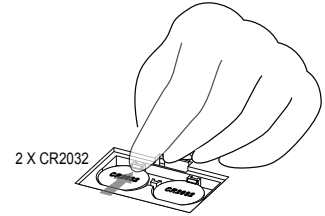*Please remove the insulating strip before use.

ADDITIONAL INFORMATION:

- 1.) Make sure your scale is at room temperature, and on a stable horizontal surface free of vibration.
- 2.) This scale is not waterproof. Do not expose to moisture.
- 3.) As this scale is a precision instrument, it must be treated with extreme care. Avoid rough treatment, shaking and vibration.
- 4.) Do not disassemble this scale or the warranty will be void.
- 5.) The display value may be affected by electromagnetic disturbances such as operating near a radio. If this occurs, move the device that is causing the disturbance and turn the scale off and restart.

PRODUCT GUARANTEE

Your Escali scale is guaranteed to be free of defects in materials and workmanship for a lifetime. For full details of this guarantee please visit our website at:

see www.escali.com for full details
© Copyright Escali LLC, 2006

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vea www.escali.com para detalles completos © Escali LLC, 2009