Animal Sounds Digital Thermometer


Item 15-901-000, 15-902-000, 15-903-000

Thank you for purchasing a HealthSmart Kids Digital Thermometer. Please read the following detailed instructions to achieve the most accurate temperatures and safe operation. With proper care and use, your thermometer will provide you with reliable readings.

What is a "Normal" Temperature?

Although the generally accepted "normal" temperature reading is 98.6°F (37.0°C), temperature readings can vary from 97°F (36.1°C) to 99°F (37.2°C) and still may be considered "normal". Variations in temperature can be attributed to activities such as exercise, smoking, eating and drinking. Even time of day may influence your temperature. For example, your temperature is lower in the morning than in the afternoon. Other variations may be due to the location of the temperature reading. While the oral temperature follows the guidelines previously described, a rectal temperature is generally 1.0°F (0.5°C) higher. Conversely, an axillary (under the arm) temperature will be 1.0°F (0.5°C) lower.

Features of your Digital Thermometer

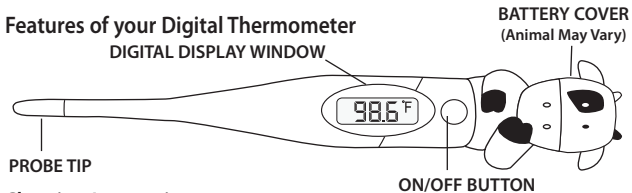Cleaning Instructions

Before and after each use, clean the thermometer with a soft cloth and 70% isopropyl alcohol diluted with water, or wash with COLD soapy water. Do not submerge the thermometer in any liquid. Wipe the thermometer with a dry cloth and allow to air dry.

1. The unit must never be sterilized by boiling, gas or steam autoclaves.
2. The use of any other solution or cleaning method beyond the above mentioned may result in reliability problems with the thermometer and void your warranty.

Switching Dual Scale

NOTE: To change between Celsius and Fahrenheit readings, hold the POWER button for 2 seconds until the "Lo" with an "m" displays and changes between °F and °C. Press the POWER button again to turn the thermometer off.

How to Take Your Temperature

1. Clean the thermometer according to the instructions provided.
2. Apply a disposable probe cover to the clean, dry thermometer.
3. Determine whether the measurement site will be oral, rectal or underarm.

NOTE: Measurement time may vary based on the selected site.

Once the desired temperature is selected, press the POWER button (the unit will beep) to activate the thermometer. This is the LCD display test, which insures all parts of the display are functioning properly.

The last measured temperature with an "m" will appear. If there is no previous measurement, "Lo" °F(°C) will be displayed.

98.6°F(37.0°C) will briefly appear

Finally, a "Lo" with a flashing °F(°C) will appear in the display to indicate that the thermometer is ready for measurement.

NOTE: If the room temperature is above 89.6°F(32.0°C), the room temperature will display instead of the Lo.

Taking an Oral Temperature


Approximate measurement time: 30 seconds

1. Using the diagram on the right, place the probe tip under the tongue as near as possible to the heat pocket. The diagram also shows the temperature variables by not placing the tip in the heat pocket area.

NOTE: Drinking hot or cold fluids, exercising, smoking or performing other activities may raise or lower your temperature. Therefore, it is important to relax for approximately 5 minutes with your mouth closed prior to taking a reading.

2. The current temperature is shown and the °F (°C) will continue to flash during the reading.
3. When the peak temperature has been reached, the °F(°C) will stop flashing and the thermometer will sound four long barks, moos or quacks if the temperature is below 100°F(37.8°C). If the temperature is above 100°F(37.8°C) the unit will sound two sets of three short barks, moos or quacks.
4. To turn the unit off, press the ON/OFF button. The unit will automatically shut off after approximately ten minutes.

NOTE: The mouth must remain closed during the measurement for the most accurate reading. Opening the mouth could result in an extended measurement time and may affect the reading.

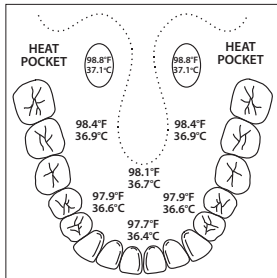Taking an Axillary Temperature (under the arm)

Approximate measurement time: 40 seconds

1. Wipe the underarm with a dry towel.
2. Place the probe tip under the arm so the tip is touching the skin with the thermometer perpendicular to the body.
3. Position the arm across the chest so the probe tip is well covered by the arm. This also insures the probe is not affected by the room's air.
4. When the peak temperature has been reached, the °F(°C) will stop flashing and the thermometer will sound four long barks, moos or quacks if the temperature is below 100°F(37.8°C). If the temperature is above 100°F(37.8°C) the unit will sound two sets of three short barks, moos or quacks.
5. To turn the unit off, press the ON/OFF button. The unit will automatically shut off after approximately ten minutes.

NOTE: Axillary temperatures are generally 1.0°F (0.5°C) lower than oral readings.

Taking a Rectal Temperature

Approximate measurement time: 20 seconds

NOTE: Rectal temperatures are an appropriate and reliable method for infants and small children. Otherwise, this method should only be used when it is impossible or impractical to take an oral or axillary temperature.

1. After applying the disposable probe cover, apply a water-soluble lubricant to the tip.
2. Gently insert the probe a MAXIMUM of 1/2" into the rectum.
3. When the peak temperature has been reached, the °F(°C) will stop flashing and the thermometer will sound four long barks, moos or quacks if the temperature is below 100°F(37.8°C). If the temperature is above 100°F(37.8°C) the unit will sound two sets of three short barks, moos or quacks.
4. To turn the unit off, press the ON/OFF button. The unit will automatically shut off after approximately ten minutes.

NOTE: Rectal temperatures are generally 1.0°F (0.5°C) higher than oral readings.

Using Disposable Probe Covers

It is recommended to use the enclosed disposable probe covers to prevent cross contamination. Additional probe covers can be purchased from your retailer or with the order form provided.

General Care

1. The thermometer should only be used under the supervision of an adult.
2. Do not walk, run or talk during temperature taking.
3. Clean the thermometer before and after each use.
4. Store the unit in the protective case when not in use.
5. Do not bite the probe.
6. Do not store the unit where it will be exposed to direct sunlight, dust or humidity. Avoid extreme temperatures.
7. Dropping or subjecting your thermometer to strong shocks should be avoided.
8. Do not attempt to disassemble the unit, except to replace the battery.
9. If measurement is not attainable, the unit must be replaced.

NOTE: Performance of the device may be degraded if: operated outside of stated temperature and humidity range; stored outside of stated temperature and humidity range; thermometer undergoes mechanical shock (drop); patient temperature is below ambient temperature.



Changing the Battery

The battery in your digital thermometer needs to be replaced when the "■" symbol appears. Replace the battery with two LR1130 or equivalent.

To replace the battery,

1. Remove screw and slide the battery cap away from the unit.
2. Gently pull the battery tray out until it stops. Using a non-metal pointed tool, carefully remove the batteries and insert the new batteries with the positive (+) side facing up.
3. Carefully slide the battery cap back into place, making sure that the cap fits snugly to protect the battery compartment from moisture and replace screw.

NOTE: Please properly dispose of the batteries away from children and heat.

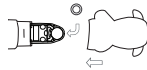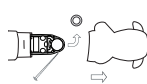Specifications

Item 15-901-000, 15-902-000, 15-903-000

DISPLAY:	Liquid Crystal Display
DISPLAY RESOLUTION:	4 Digits (1/10TH of a degree increments)
MEASUREMENT RANGE:	90.0°F to 109.9°F (32°C - 42.9°C)
ACCURACY:	±0.2°F (98.0°F - 102.0°F) ±0.1°C (37.0°C - 39.0°C)
AVERAGE MEASUREMENT TIME:	Oral 30 seconds. Rectal 20 seconds. Axillary 40 seconds
TO NE:	Thermometer barks, moos or quacks when activated and deactivated. A series of barks, moos or quacks signal when peak temperature is reached.
MEMORY:	Displays last measured temperature
BATTERY:	Two, LR1130 or equivalent
BATTERY LIFE:	Approximately 1,500 measurements
AUTOMATIC SHUT-OFF:	After approximately 10 minutes
AMBIENT OPERATING TEMPERATURE:	50.0°F to 113.0°F (10°C to 45°C)
HUMIDITY:	15% to 95% RH
STORAGE TEMPERATURE:	-4.0°F to 122.0°F (-20°C to 50°C)
HUMIDITY:	15% to 95% RH
DIMENSIONS:	Dog: 5.9" x .98" x 1.02" Cow: 5.7" x 1.14" x .83" Duck: 5.8" x .94" x .91"
WEIGHT:	Dog: 0.85 ounces (24 grams) Cow: 0.81 ounces (23 grams) Duck: 0.77 ounces (22 grams)

Specifications are subject to change without notice.

#91-032-901 01/11

©2011 Briggs Medical Service Company

LIMITED FIVE-YEAR WARRANTY

The warrantor guarantees that this Thermometer will be free from manufacturing defects under normal use for five years from the original purchase date. This warranty covers only normal use, and does not apply to use in any clinical or commercial applications. This warranty does not cover batteries that may be provided with or used with the Thermometer. This warranty is voided if the Thermometer product is misused or abused in any manner.

If the Thermometer fails to operate during the first five years from the time originally purchased, please call our Customer Care Help Line at 800-622-4714 and we will try to resolve the issue as soon as possible. If the warrantor determines the unit failed to operate due to a manufacturing defect, the unit will be replaced at the option of the warrantor. Replacement is the sole remedy under this limited warranty. This warranty gives you specific legal rights, which vary from state to state. As a condition of this warranty, the enclosed warranty registration card must be completed and sent to us within 10 days of the purchase date or can be filled out online at www.mabisdmi.com.

This warranty constitutes the warrantor's only responsibility and obligation to replace materials or components. We make no other express or implied warranties, arising by operation of law or otherwise, or any warranty of merchantability or fitness for a particular use or purpose whether or not the use or purpose has been disclosed to the warrantor in specifications, drawings or otherwise, and whether or not the warrantor's products are specifically designed and/or manufactured by the warrantor for the buyer's use or purposes, except for the limited warranty stated above. The warrantor will not be responsible for any indirect, incidental, special, consequential, or punitive damages or other loss, including, but not limited to, damage to or loss of other property or equipment and personal injuries, whether to purchaser or others. The warrantor shall in no event be liable to the purchaser for any amount in excess of the cost of replacement of the unit.

CUT OUT AND MAIL TO:
HEALTHSMART
1931 NORMAN DRIVE SOUTH
WAUKEGAN, IL 60085

THERMOMETER WARRANTY REGISTRATION

(MUST BE COMPLETED WITHIN 10 DAYS OF PURCHASE)

NAME: _____

ADDRESS: _____

CITY: _____ STATE: _____ ZIP: _____

PURCHASED FROM: _____

CITY: _____ STATE: _____ ZIP: _____

DATE PURCHASED: _____

MODEL PURCHASED: 15-901-000 15-902-000 15-903-000

USER INFORMATION (OPTIONAL)

MALE FEMALE AGE: _____

AGE OF YOUNGEST FAMILY MEMBER

UNDER 6 MONTHS 6 MONTHS - 1 YEAR
 1 - 3 YEARS 3 - 5 YEARS OLDER

PROBE COVER ORDER FORM

CUT OUT AND MAIL TO:

HEALTHSMART
1931 NORMAN DRIVE SOUTH
WAUKEGAN, IL 60085

Please send _____ packages of disposable probe covers #15-618-000 (100 probes per package). Enclosed is \$5.00 per package plus \$3.00 for shipping and handling for each order.

PAYMENT METHOD (do not send cash): Check Money Order

Please ship to:

NAME: _____

ADDRESS: _____

CITY: _____ STATE: _____ ZIP: _____

PHONE: (_____) _____

NOTE: This order form is for shipments within the United States and Canada.

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