

## LogicMark

## Quick Start Guide




STEP 1: Plug phone cord and power cord (AC Adapter) into the back of the base. Give power cord a half turn to the right to lock in place.

STEP 2: Plug phone cord and power cord (AC Adapter) into wall jacks.



If you have any questions regarding this product, please call: 1-800-519-2419 1-703-934-7934

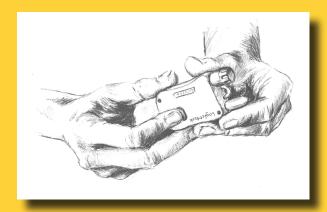

STEP 3: Run a full system check to make sure unit is working. Press the grey button on the back of the pendant and hold for <u>four full seconds</u>. Unit will beep, and after a few seconds you should hear a dial tone and the announcement "All systems are ok."