



SUNRISE MEDICAL improving people's lives User Name:dealerdemo Sold To:684741

My C2 Address Book Inquiry Orders/Quotes Inquiries RMAs Reports

Start New Order Start New Quotes Edit/Release Order/Quotes Begin New Sample

Shopping Cart

Option Selection

M2100 -- (1/EA)

Please enter information and select Add or Cancel to continue.
Note: Additional days may be applied, please reference the schedule date on your order

CANCEL >> **ADD >>**

<p>*JYW-SIZE WIDTH-FG9-SIZE WIDTH</p> <p>Please Select</p>	<p>*JYW-SIZE LENGTH-FG2-SIZE LENGTH</p> <p>Please Select</p> <p>JYW-IL-OPT01 - IL1 -INCREASE LENGTH 1" JYW-IL-OPT02 - IL2 -INCREASE LENGTH 2" JYW-IL-OPT03 - IL3 -INCREASE LENGTH 3" JYW-IL-OPT04 - IL4 -INCREASE LENGTH 4" JYW-IL-OPT05 - IL5 -INCREASE LENGTH 5" JYW-RL-OPT01 - RL1 -REDUCE LENGTH 1" JYW-RL-OPT02 - RL2 -REDUCE LENGTH 2" JYW-RL-OPT03 - RL3 -REDUCE LENGTH 3" JYW-RL-OPT04 - RL4 -REDUCE LENGTH 4" JYW-RL-OPT05 - RL5 -REDUCE LENGTH 5"</p>
<p>*JYW-SIZE SLOPE-FG2-SIZE SLOPE</p> <p>Please Select</p>	<p>*JYW-LAT PELVIC-FG2-LATERAL PELVIC POSITION</p> <p>Please Select</p>
<p>*JYW-FLUID OB LT-FG1-FLUID PAD OBLIQUITY LEFT</p> <p>Please Select</p>	<p>*JYW-COVER EMBR-FG1-COVER EMBROIDERY</p> <p>Please Select</p>
<p>*JYW-BASE MISC-FG3-BASE MISCELLANEOUS</p> <p>Please Select</p>	
<p>*JYW-LEG TR RT-FG2-LEG TROUGH CUT RIGHT</p> <p>Please Select</p>	
<p>*JYW-MEDIAL THIGH-FG2-MEDIAL THIGH POSITION</p> <p>Please Select</p>	
<p>*JYW-COVER MAT-FG2-COVER MATERIALS</p> <p>Please Select</p>	
<p>*JYW-SIDE COLOR-FG1-COVER SIDEBANDS COLOR</p> <p>Please Select</p>	

Select JAY cushions and backs are now easier than ever to order with JAY Your Way modifications using our online ordering process! Follow these simple steps to ensure you receive a unique JAY product that is manufactured to your exact specifications. For an updated list of cushions available to order online with modifications, please visit the JAY Your Way product page on www.sunrisemedical.com.

- Step 1:** Consult the JAY Your Way Catalog and JAY Your Way Price List to select the modifications you need.
- Step 2:** Visit www.sunrisemedical.com
- Step 3:** Log on to your My Sunrise account.*
- Step 4:** Begin an order or quote.
- Step 5:** Enter an "M", followed by the item number of the cushion or back you would like to modify.
- Step 6:** Select from drop-down menus the modifications you require.


*To set up a My Sunrise account, visit www.sunrisemedical.com, click "Sign up for an Online Account," and follow the prompts.



JAY Your Way is a modification program designed to individualize and customize JAY products. From clinical improvements to colored covers, JAY Your Way can be used to create cushions and backs manufactured to users' unique individual needs. This catalog includes information on customizing JAY bases, fluid inserts, covers and backrests, and using our simple online ordering process.

Cushion Base Modifications


Size Modifications

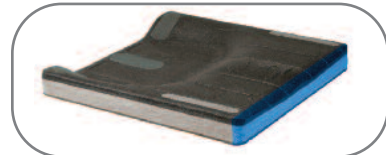

Reducing Size

Mod Detail: Removing foam from the length or width of the cushion base and sewing the cover to fit

Clinical Application: To match the size of the client and/or wheelchair

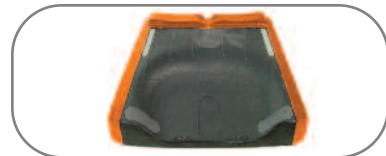

Increasing Size

Mod Detail: Adding dummy foam to the length or width of the cushion base and sewing the cover to fit


Clinical Application: To match the size of the client and/or wheelchair



Increasing Height

Mod Detail: Adding foam underneath the cushion base to increase height/thickness

Clinical Application: To accommodate wheelchair fitting and/or environment needs e.g transfers or ground clearance

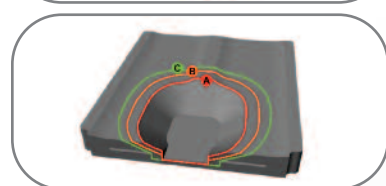

Well Sizing

Mod Detail: Increasing or decreasing the standard, factory set size of the Pelvic Loading Area (PLA)

Clinical Application: To allow for anatomically optimized loading areas when considering proper immersion, pressure reduction and lateral stability

Base and Positioning Modifications



Notch for Rails

Mod Detail: Removing foam along the bottom, lateral edges of the foam base


Clinical Application: To lower the cushion between seat rails to achieve the lowest possible seat-to-floor height



Notch for Posts

Mod Detail: Notching the back corners of the foam base


Clinical Application: To slide cushion between back posts when loading surface is posterior to back posts, to lock in cushion in the case of frequent transfers and to build in some depth growth opportunities

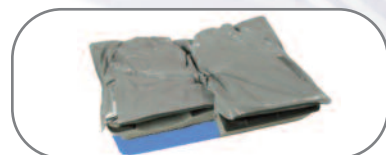

Leg Length Cut (Right/Left)

Mod Detail: Decreasing the length of the foam base and fluid pad (if applicable) on the left or right side and sewing the cover to fit

Clinical Application: To accommodate varying femoral lengths

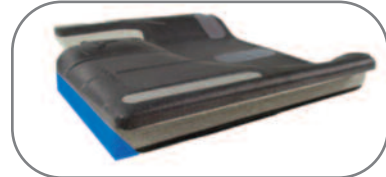

Contracture Cut

Mod Detail: Undercutting the front of the foam base at an angle


Clinical Application: To clear the calves in cases of tight hamstrings, large calves, or close front thigh to lower leg angles, to accommodate chair configurations, and to accommodate foot propulsion



Scrotal Relief

Mod Detail: Scooping out the front and center of the Pelvic Loading Area (PLA)


Clinical Application: To accommodate organ tissue, reduce discomfort, and increase sitting tolerance

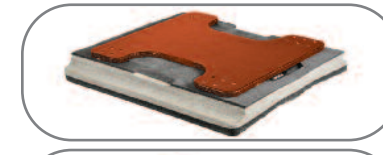

Custom Drop Seat

Mod Detail: Fitting a drop seat to the cushion base

Clinical Application: To minimize seat-to-floor height for foot propelling and transfers

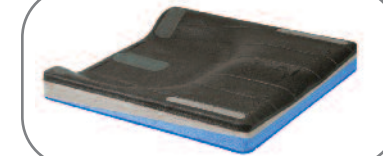

Angle Cut Posterior/Anterior Slope

Mod Detail: Removing foam at an angle from underneath the cushion base

Clinical Application: **Anterior Slope:** To accommodate hip flexion limitations and assist with foot propelling and transfers. **Posterior Slope:** To close seat to back angle (e.g. for postural stability)



Wedge Glue Posterior/Anterior Slope

Mod Detail: Adding a foam wedge underneath the cushion base

Clinical Application: **Anterior Slope:** To accommodate hip flexion limitations and assist with foot propelling and transfers. **Posterior Slope:** To close seat to back angle (e.g. for postural stability)



Medial Support Cut

Mod Detail: Removing the pre-contoured medial support

Clinical Application: To improve positioning when femoral alignment is not neutral (e.g. in windswept postures) or when transfers are inhibited by the medial thigh support



Femoral Relief

Mod Detail: Shaving away the front of PLA to transition from PLA to femoral support area


Clinical Application: To reduce pressure at front of PLA area to accommodate tissue and increase sitting tolerance for users with sensation



Leg Trough Cut

Mod Detail: Scooping away 1" of base material from the front of the PLA forward to the front edge of the femoral support area

Clinical Application: To improve medial and lateral positioning and support of the femurs and to accommodate hip limitations by opening hip angles

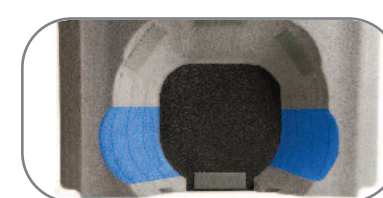


Trochanteric Relief

Mod Detail: Removing material at the lateral borders of the PLA, beneath the trochanters


Clinical Application: To reduce pressure at or surrounding the trochanters and evenly distribute load across the whole sitting surface

Positioning Components

Positioning pieces, such as adductors, abductors, and hip guides, can be added to several foam cushion bases to improve sitting posture and positioning. The below positioning components are shown on the JAY J3 Cushion base.

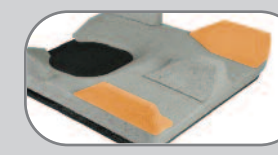

Lateral Thigh

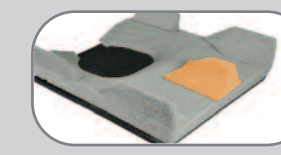

Medial Thigh



Lateral Pelvic

Fluid Insert Modifications



Underfill

Mod Detail: Reducing factory set fluid volume in the JAY Flow fluid pad

Clinical Application: To increase immersion of the pelvis into the fluid to achieve appropriate skin protection and/or lateral stability when bottoming out is not a concern



Overfill

Mod Detail: Increasing factory set fluid volume in the JAY Flow fluid pad


Clinical Application: To redistribute load and reduce pressure on the bony prominences when a client experiences "bottoming out" or discomfort on a standard JAY Flow fluid pad



Accommodating Obliquities

Mod Detail: Applying an overfill to either the right or left side of the JAY Flow fluid pad

Clinical Application: To optimize immersion (underfill), prevent bottoming out (overfill), or substitute for bone or tissue asymmetry (overfill) on one side of the client in the case of a flexible obliquity

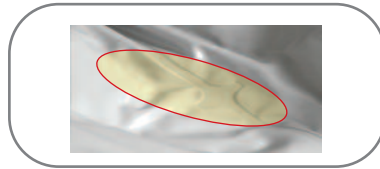

Sacral or Scrotal Seal

Mod Detail: Adding a heat seal in the sacral or scrotal area of the JAY Flow fluid pad

Clinical Application: To control the migration of fluid in order to reduce pressure in areas that are at-risk for skin breakdown, pain, or discomfort

Cover Modifications

Each JAY cushion and back comes with a carefully designed cover manufactured from top-of-the-line materials that enhance their clinical effectiveness and comfort. JAY Your Way can be used to add some color or better fit users' unique lifestyles.



Sealed Inner Cover

Mod Detail: The water-resistant, stretchable Sealed Inner Cover reduces surface tension and protects from dirt, moisture and odor build-up.

Clinical Application: To reduce surface tension, to protect foam base from dirt and moisture build-up, and to prevent odor



Custom Outer Covers

Specialty outer covers are available on select JAY cushions and backs. Options include:

- Incontinent Shower Cap Cover
- Incontinent Zipper Cover
- JAY Stretch Zipper Cover
- X-static® Cover with Breathable Spacer Fabric
- Dartex™ Cover with Aquaguard Zipper
- Reverse Dartex™ Cover with Aquaguard Zipper



Embroidery

Adding embroidery to JAY cushions and backs is a great way to personalize or keep track of products in larger facilities. Available in red, blue, pink, purple, toxic green, forest green, orange, gray, white or black Arial font



Colored Outer Cover with Zipper

Select JAY cushions are available with colored outer covers or colored semi-stretch sidebands for added flair! Available color options are red, blue pink, purple and toxic green.

Backrest Modifications



Military Service Patches

To show past service in the military or support for the armed forces, Military Service Patches can be added to the JAY J3 and GO backs. Military Service Patches are also available on Quickie wheelchairs' tension adjustable upholstery!



JAY J3 Back Cover


A sleek black or rugged desert camouflage cover can be chosen at no charge to replace the JAY J3's standard gray cover.



Spinal Fluid Pad

Mod Detail: Adding a JAY Flow insert over spinal area of foam back

Clinical Application: To reduce pressure along the spine, which is often at high risk of skin breakdown due to its bony structure

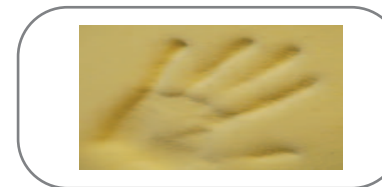

Visco Foam Insert

Mod Detail: Adding a Visco memory foam insert with a cover sewn to match

Clinical Application: To add a comfort layer over the standard supportive foam for increased sitting tolerance and pressure reduction



Painted Back Shell

Select JAY backrests are available with painted back shells to match or complement the wheelchair! Paint colors include glossy Candy Red, Candy Blue, Forest Green, Silver Texture, Pink, Black, and matte Black Slate.

